

## Product Specification

# Prep and Pack Ergonomic Workstations

**Designed with clinician comfort and safety in mind, our prep and pack ergonomic workstations are ideal for a variety of instrument processing tasks.**

### Product Description

Belimed's line of Prep & Pack workstations are designed to promote efficient, ergonomic processing and packaging of instruments and utensils. Workstations are highly configurable, available in various sizes, with different work surface materials, and with a wide array of storage, lighting, and technology support options. Frames are also available for each size, in fixed, manual adjust, and power adjustable configurations.

### Application

For use in healthcare facilities for the inspection, organization, and wrapping of re-usable instruments and utensils prior to sterilization.

### Standard Features

#### Construction/Design

Belimed Prep & Pack workstations are constructed of high-quality materials and finishes and include the following standard features:

- Type 304 stainless steel standard tops with seam welded corners and vibration-damping solid core
- Heavy-gauge steel frames with durable powder-coat paint surfaces
- Cantilevered frame design (powered & manual) for maximum knee well access
- Up to 1,200 lb. lift capacity on electric adjust units
- Slotted uprights for tool-free overhead storage reconfiguration
- Widest range of accessories

### Benefits

#### Optimized Ergonomics

Manual height adjustment (ideal for the single shift/operator station), and electric powered height adjustment (best choice for multiple operator/shift applications) are available. Programmable controls that store user height settings are a standard option. Fixed-height is also available as an entry-level option. All workstations and accessories are designed to adjust to fit the needs of 5th through 95th percentile users, having been built to the ANSI/HFS 100-1988 standard.

#### Ideal for Multiple Users

Power adjust models feature smooth action with simple easy to use under-counter up/down controls. Each user may select the proper working height to fit their stature, thereby ensuring hours of fatigue-free operation. A height-adjustable workstation allows the users to adjust to an optimized, neutral position and reduce or eliminate the ergonomically risky and non-value adding body movements while working.

#### Massive Stability

All frames are fully welded of heavy-gauge steel and include back rails to prevent frame distortion and racking under load. Solid core stainless tops reduce vibration and sagging common to standard welded-gusset reinforced-top designs.

#### Easy Reconfiguration

Tool-free overhead storage system connections allow shelves and bins to be quickly and easily reconfigured, to meet changing department or user needs.

#### Serviceability

Motors and controls for power adjustable units drop out for quick and easy change-out with minimal downtime.



### Options

#### Surface Material:

- Stainless Steel
- T-Mold Plastic/Cleanroom Laminate

#### Frames:

- Power
- Manual
- Fixed

#### Workstation Widths:

- 48"
- 60"
- 72"

*More sizes available on request*

#### Work Surface Depths:

- 30"
- 36"

*More sizes available on request*

#### Wide Range of Storage & Organizational Options:

Shelves, Bin racks, Drawers, Wrap rods, Rails, and Baskets are all available to keep the supplies you need close at hand.

#### NEW OPTION – Social Distancing Dividers:

These can be used to quickly create separate workspaces. These can form a single separation wall, an L-shaped enclosure, a U-shaped enclosure, and many more configurations. The polycarbonate panel allows users to see through and it can be easily cleaned, supporting your sanitation protocol.

#### Full Technology Support:

PC's, printers, and scanners slot easily into Belimed's workstation, with mounting options for flat screens, CPU's, printers, and keyboards, and comprehensive power outlet and cord management support. Auxiliary carts printer/PC carts are also available to help maximize countertop and under counter space.



For more information, visit [belimed.com](http://belimed.com)

# Prep and Pack Ergonomic Workstations

### Accessories

#### Overhead Storage:

- 8" x 30" shelf, 12" x 30" shelf, 15" x 30" shelf
- 4" x 30" bin rail
- Spool holder rod
- Small roll shelf
- Peel pack box dispenser (for bin rail)
- Wrap paper backframe (4 rods above work surface)

#### Under-counter Storage:

- Single 6" drawer, Single 12" drawer
- Dual-stacked 6" drawers
- Dual slide-out wire baskets

#### Power Management:

- 8 outlet power bar
- Rear wire trough
- Horizontal power bar work surface mount

#### PC Support:

- Articulating keyboard arm with mouse tray (fixed or adjustable height)
- Articulating flat-panel display arm (fixed or adjustable height)
- Fixed shelf mount keyboard tray
- Pull-out/Under-surface keyboard tray
- CPU undermount bracket

#### Overhead Lighting:

- 48" Dual-tube
- 36" Dual-tube

#### Task Lighting:

- Fluorescent bar light
- Fluorescent light ring magnifier
- LED light ring magnifier

#### Organizational Aids:


- Magnetic marker whiteboard
- Pull-out work surface
- Pegboard

#### Miscellaneous:

- Adjustable footrest
- Counter backplash
- Locking & non-locking casters
- Leveling casters
- Printer carts
- Wrap paper rack
- Tool balancer trolley




8 Outlet Power Bar

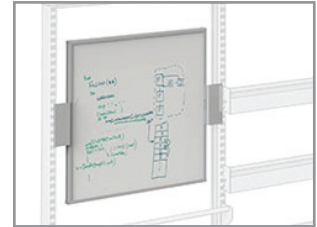


Magnetic Whiteboard



Flat Panel Display Arm

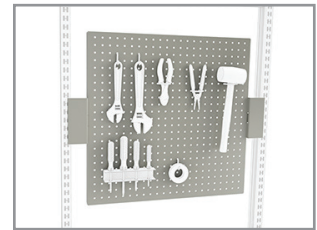


Pegboard

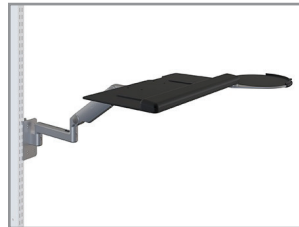

Keyboard Arm with Mouse Tray



Free Standing Footrest

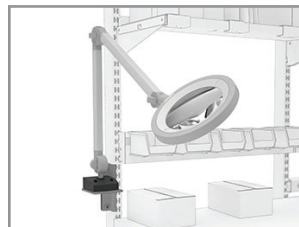

Light Ring Magnifier



Tool Balancer Trolley

Distributed by Belimed for BOSTONtec

REV5.20.2021\_Q0

Engineers of Confidence.