



EQUIPMENT ANCHORAGE & SEISMIC ENGINEERING

2801 Connery Way, Suite B

Missoula, MT 59808

Phn: (406) 541-EASE (3273) Fax: (406) 541-3274

Sheet 1 of 8

Office of Statewide Health Planning and Development

ANCHORAGE PRE-APPROVAL

OPA-2750-10

THIS PRE- APPROVAL CONFORMS TO THE 2010 CALIFORNIA BUILDING CODE

Equipment Manufacturer: Belimed Inc

Equipment Type: CS750 Cart Washer/Tank

GENERAL NOTES

1. EXPANSION ANCHORS:


- (a) ATTACHMENT IS TO BE MADE WITH THE ANCHORS LISTED BELOW AND INSTALLED AS DESCRIBED IN THE CORRESPONDING ICC REPORT.

Anchor Diameter	Concrete Type	Min. f'c (psi)	Anchor Type	ICC Report No.	Min. Embed.	Min. Spacing	Min. Edge Dist.	Min. Conc. Thickness	Installation. Torque	Test Loads
3/4"	Hardrock	3000	Hilti Kwik Bolt TZ	ESR-1917	3-3/4"	18"	33"	6"	110 Ft-Lbs	Direct Pull Tension - 4120 lbs
1/2"	Hardrock	2500	Hilti Kwik Bolt TZ	ESR-1917	2"	12"	18"	4"	40 Ft-Lbs	Direct Pull Tension - 1465 lbs

- (b) THIS PRE-APPROVAL ALLOWS FOR UP TO A MAXIMUM OF 2 ADJACENT CONCRETE SLAB EDGES, 18" OR 33" (SEE SCHEDULE) AWAY MINIMUM (i.e. - CORNER). SEE ADJACENT DETAIL FOR ADDITIONAL MINIMUM ALLOWABLE CONCRETE EDGE DISTANCES.

2. TESTING OF EXPANSION ANCHORS PER 2010 CBC, 1916A.7 : TENSION TESTING SHALL BE DONE IN THE PRESENCE OF THE SPECIAL INSPECTOR AND A REPORT OF THE TEST RESULTS SHALL BE SUBMITTED TO OSHPD

- (a) AFTER AT LEAST 24 HOURS HAVE ELAPSED SINCE INSTALLATION, DIRECT PULL TENSION TEST AT LEAST 50% OF THE ANCHORS.
- (b) ACCEPTANCE CRITERIA: THE ANCHOR SHOULD HAVE NO OBSERVABLE MOVEMENT AT THE TEST LOAD. A PRACTICAL WAY TO DETERMINE OBSERVABLE MOVEMENT IS THAT THE WASHER BECOMES LOOSE.
- (c) IF ANY ANCHOR FAILS, TEST ALL ANCHORS.

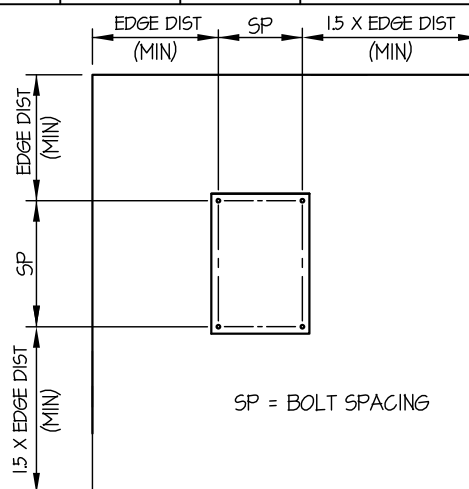

TYPICAL CONCRETE EDGE DETAIL



A P P R O V E D	
Fixed Equipment Anchorage Office of Statewide Health Planning and Development	
OPA-2750-10	
Pre-approval Program Manager: Anthony R. Pike (916) 440-8470	
P E N D I N G	
Reviewed By: Anthony R. Pike	9/22/11



EQUIPMENT ANCHORAGE & SEISMIC ENGINEERING

2801 Connery Way, Suite B

Missoula, MT 59808

Phn: (406) 541-EASE (3273) Fax: (406) 541-3274

Sheet 2 of 8

Office of Statewide Health Planning and Development
ANCHORAGE PRE-APPROVAL

OPA-2750-10

THIS PRE- APPROVAL CONFORMS TO THE 2010 CALIFORNIA BUILDING CODE

Equipment Manufacturer: Belimed Inc

Equipment Type: CS750 Cart Washer/Tank

GENERAL NOTES (CONTINUED)

3. FORCES PER ASCE 7-05 SECTION 13.3.1, EQUATIONS 13.3-1, 13.3-2 & 13.3-3, WHERE $S_{Ds} = 2.0$, $a_p = 2.5$, $I_p = 1.5$ & $R_p = 6.0$, $z/h = 0.0$ AT CONCRETE SLAB & $z/h = 1.0$ AT CONCRETE ON METAL DECK.
4. THIS PRE-APPROVAL CONFORMS TO THE 2010 CALIFORNIA BUILDING CODE.
5. THIS PRE-APPROVAL COVER ONLY THE ANCHORAGE OF THE EQUIPMENT TO THE HOSPITAL BUILDING'S STRUCTURE.
6. ALL ANCHOR FORCES SHOWN ON THE DRAWINGS ARE FACTORED LOADS THAT SHALL BE USED FOR STRENGTH DESIGN.

RESPONSIBILITIES OF THE STRUCTURAL ENGINEER OF RECORD OF THE BUILDING

7. VERIFY THAT PROJECT SPECIFIC VALUES OF S_{Ds} & z/h RESULT IN SEISMIC FORCES (E_h , E_v) THAT DO NOT EXCEED THE VALUES ON THE DETAILS.
8. VERIFY THAT THE CONCRETE SLAB TO WHICH THE EQUIPMENT IS ANCHORED MEETS THE REQUIREMENTS OF THE APPLICABLE ICC ESR
9. VERIFY THAT THE ANCHORS ARE AN ADEQUATE DISTANCE FROM ANY SLAB EDGES OR OPENINGS (SEE TYPICAL DETAIL ON SHEET 1).
10. VERIFY THAT ALL NEW OR EXISTING ANCHORS ARE AN ADEQUATE DISTANCE FROM THE ANCHORS SHOWN IN THIS PRE-APPROVAL. SEOR SHALL VERIFY THAT THERE IS NO ADVERSE INTERACTION WHERE OTHER ANCHORS ARE WITHIN 18" OR $6h_{ef}$ FROM THIS UNIT'S ANCHORS
11. DESIGN BACKING PLATES, STUDS, ETC TO WHICH THE EQUIPMENT IS ATTACHED, AS NOTED ON THE DRAWINGS.
12. PROVIDE SUPPORTING STRUCTURE TO SUPPORT WEIGHTS AND FORCES SHOWN, IN ADDITION TO ALL OTHER LOADS. VERIFY THE ADEQUACY OF THE STRUCTURES (SUCH AS WALLS AND FLOORS) WHICH SUPPORT THE EQUIPMENT FOR THE LOADS IMPOSED ON THEM BY THE EQUIPMENT IN ADDITION TO ALL OTHER LOADS.
13. VERIFY THAT THE INSTALLATION IS IN CONFORMANCE WITH THE 2010 CBC AND WITH THE DETAILS SHOWN IN THIS PRE-APPROVAL. VERIFY THAT THE ACTUAL EQUIPMENT'S WEIGHT, CG LOCATION, ANCHOR LOCATIONS, ANCHOR DETAILS, AND THE MATERIAL AND GAGE OF THE UNIT WHERE ATTACHMENTS ARE MADE AGREE WITH THE INFORMATION SHOWN ON THE PRE-APPROVAL DOCUMENTS.



A P P R O V E D	
Fixed Equipment Anchorage Office of Statewide Health Planning and Development	
	OPA-2750-10
Pre-approval Program Manager: Anthony R. Pike (916) 440-8470	
P E N D I N G	
Reviewed By: Anthony R. Pike	9/22/11

BELIMED, INC.

CS750 CART WASHER

DES. **J. ROBERSON**

JOB NO. **11-1170**

DATE **9/22/11**


SHEET

3

OF **8** SHEETS

SEISMIC ANCHORAGE

CONCRETE SLAB

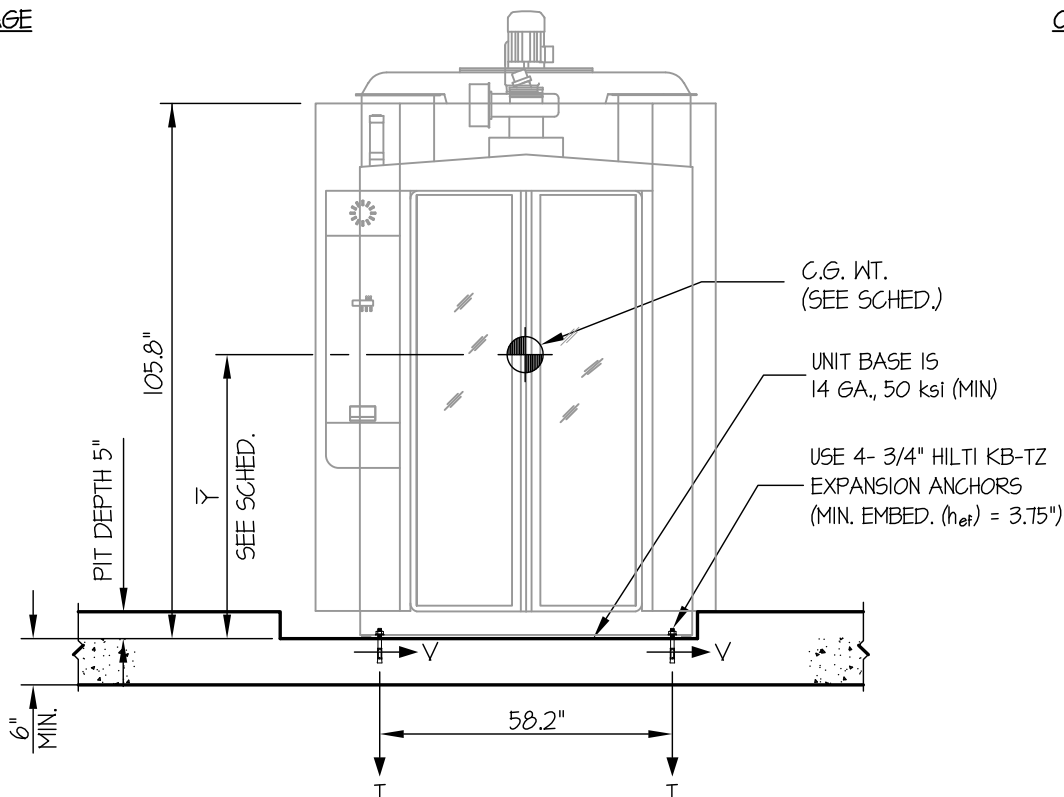

ELEVATION

SEE SHEET 2 OF 3 FOR
SCHEDULE

NOTES:

- FORCES ARE DETERMINED PER 2010 CALIFORNIA BUILDING CODE SECTION 1613A AND ASCE 7-08 SECTIONS 12 AND 13. STRENGTH DESIGN IS USED.

HORIZONTAL FORCE (E_h) = $0.90 W_p$ ($S_{DS} = 2.00$, $a_p = 2.5$, $I_p = 1.5$, $R_p = 6.0$, $z/h = 0.0$)

VERTICAL FORCE (E_v) = $0.40 W_p$

- CENTER OF GRAVITY (C.G.) WEIGHT IS A MAXIMUM. THIS CALCULATION ENCOMPASSES ALL WEIGHTS UP TO THE MAXIMUM WEIGHT SHOWN.
- ENGINEER OF RECORD FOR THE BUILDING SHALL PROVIDE SUPPORT STRUCTURE DESIGNED TO SUPPORT WEIGHTS AND FORCES SHOWN, IN ADDITION TO ALL OTHER LOADS.



APPROVED
Fixed Equipment Anchorage
Office of Statewide Health Planning and Development



OPA-2750-10

Pre-approval Program Manager:
Anthony R. Pike
(916) 440-8470

P E N D I N G

Reviewed By: Anthony R. Pike

9/22/11



EQUIPMENT ANCHORAGE & SEISMIC ENGINEERING

www.equipmentanchorage.com

BELIMED, INC.

CS750 CART WASHER

DES. J. ROBERSON

JOB NO. 11-1170

DATE 9/22/11

SHEET

4


OF

8

SHEETS

SEISMIC ANCHORAGE

CONCRETE SLAB

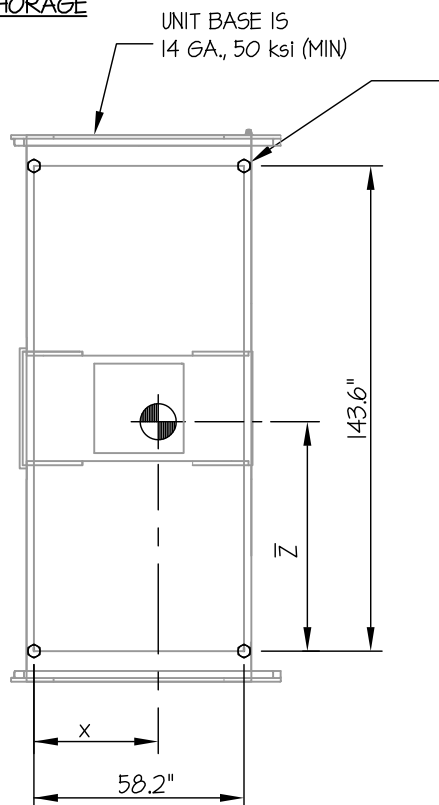

PLAN AT BASE

USE 4- 3/4" HILTI KB-TZ
EXPANSION ANCHORS
(MIN. EMBED. (h_{ef}) = 3.75")

MODEL	WEIGHT (LBS)	DIMENSIONS (INCHES)			BOLT FORCES	
		\bar{X}	\bar{Z}	\bar{Y}	T_u (LBS)	V_u (LBS)
* CS750-L	6700	33.4	63.1	53.5	2594	1730
CS750-M	5800	34.2	69.6	55.3	2188	1534
CS750-S	5400	19.3	39.8	21.3	2186	1361

* THIS MODEL REPRESENTED IN CALCULATION BELOW

 $T_u = 2594$ LB/BOLT (MAX) $V_u = 1730$ LB/BOLT (MAX)

APPROVED

Fixed Equipment Anchorage
Office of Statewide Health Planning and Development

OPA-2750-10

Pre-approval Program Manager:
Anthony R. Pike
(916) 440-8470

P E N D I N G

Reviewed By: Anthony R. Pike

9/22/11

BELIMED, INC.

CS750 CART WASHER TANK

DES. **J. ROBERSON**

JOB NO. **11-1170**

DATE **9/22/11**



SHEET

5

OF **8** SHEETS

SEISMIC ANCHORAGE

CONCRETE SLAB

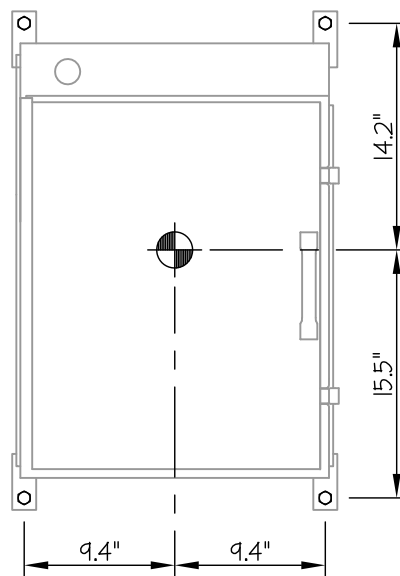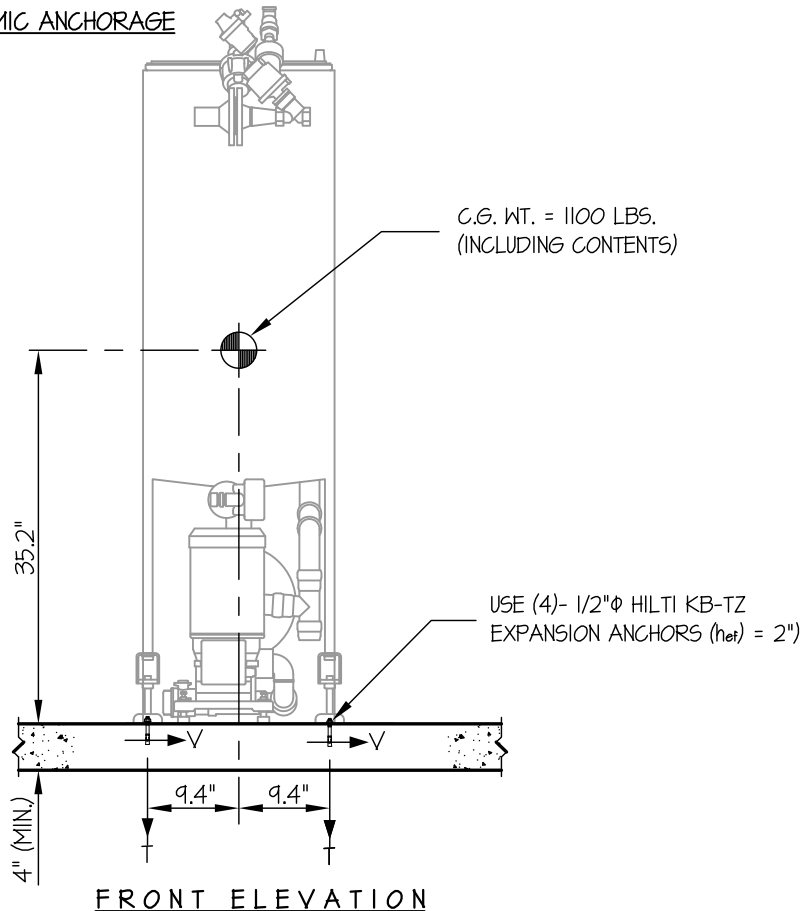

PLAN AT BASE

T_u = 1000 LB/BOLT (MAX)

V_u = 258 LB/BOLT (MAX)



A P P R O V E D

Fixed Equipment Anchorage
Office of Statewide Health Planning and Development



OPA-2750-10

Pre-approval Program Manager:
Anthony R. Pike
(916) 440-8470

P E N D I N G

Reviewed By: Anthony R. Pike

9/22/11

BELIMED, INC.

CS750 CART WASHER

DES. **J. ROBERSON**

JOB NO. **11-1170**

DATE **9/22/11**


SHEET

6

OF **8** SHEETS

SEISMIC ANCHORAGE

CONCRETE SLAB ON METAL DECK

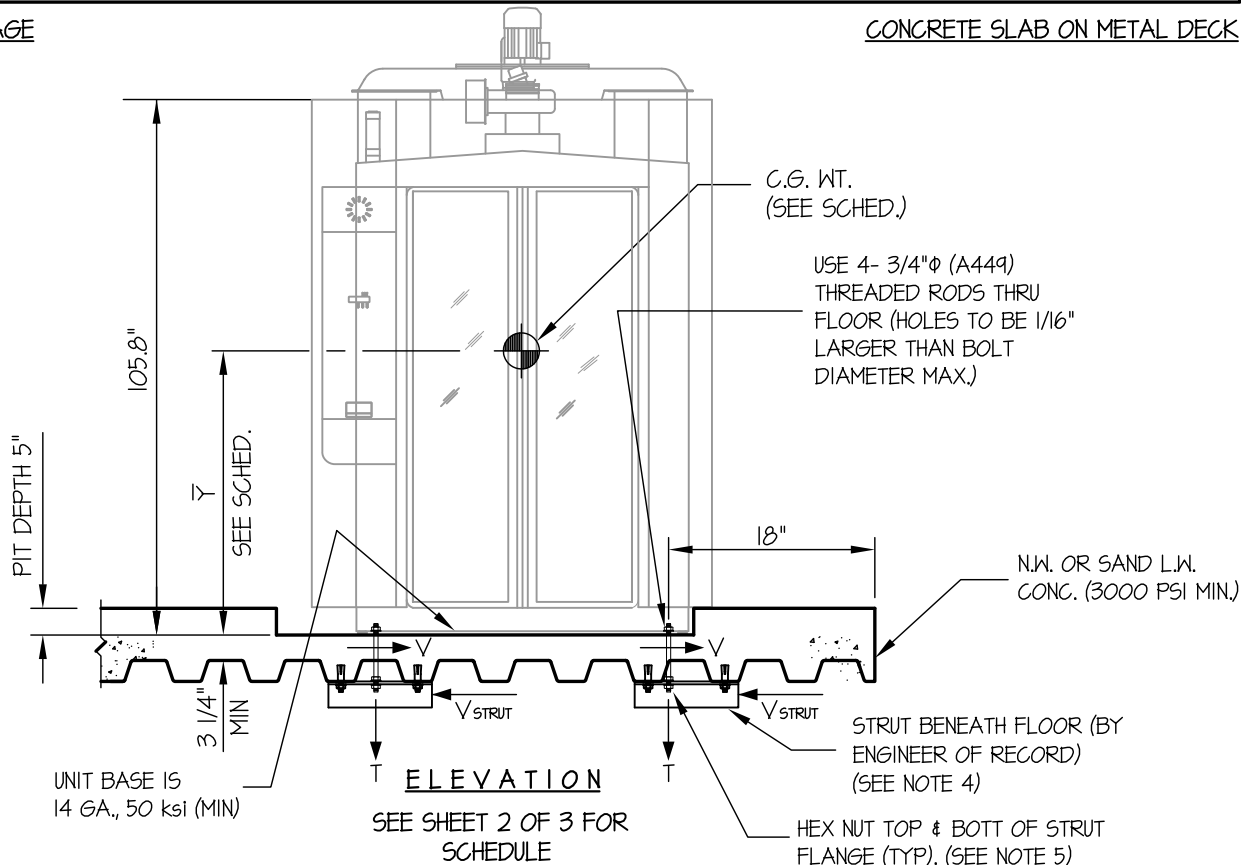

NOTES:

- FORCES ARE DETERMINED PER 2010 CALIFORNIA BUILDING CODE SECTION 1613A AND ASCE 7-05 SECTIONS 12 AND 13. STRENGTH DESIGN IS USED.
HORIZONTAL FORCE (E_h) = $1.50 W_p$ ($S_{ds} = 2.0$, $a_p = 2.5$, $I_p = 1.5$, $R_p = 6.0$, $z/h \leq 1.0$)
VERTICAL FORCE (E_v) = $0.40 W_p$
- CENTER OF GRAVITY (C.G.) WEIGHT IS A MAXIMUM. THIS CALCULATION ENCOMPASSES ALL WEIGHTS UP TO THE MAXIMUM WEIGHT SHOWN.
- ARCHITECT OR STRUCTURAL ENGINEER OF RECORD SHALL PROVIDE STRUCTURE TO SUPPORT WEIGHTS AND FORCES SHOWN.
- ARCHITECT OR STRUCTURAL ENGINEER OF RECORD SHALL DESIGN THE STRUT(S) AND ITS ATTACHMENTS TO RESIST A LOAD NOT LESS THAN V_{STRUT} IN COMBINATION WITH ALL OTHER LOADS THAT MAY BE PRESENT, WHERE $V_{STRUT} = 0.6 V_{MAX} \times$ (NO. OF ANCHORS ENGAGED BY STRUT) (MIN)
- AT CONDITIONS WHERE NUT CANNOT BE PROVIDED AT TOP SIDE OF STRUT, PROVIDE TAPPED HOLE THROUGH STRUT FLANGE.



A P P R O V E D	
Fixed Equipment Anchorage Office of Statewide Health Planning and Development	
OPA-2750-10	
Pre-approval Program Manager: Anthony R. Pike (916) 440-8470	
P E N D I N G	
Reviewed By: Anthony R. Pike	9/22/11

BELIMED, INC.

CS750 CART WASHER

DES. **J. ROBERSON**

JOB NO. **11-1170**

DATE **9/22/11**


SHEET

7

OF **8** SHEETS

SEISMIC ANCHORAGE

CONCRETE SLAB ON METAL DECK

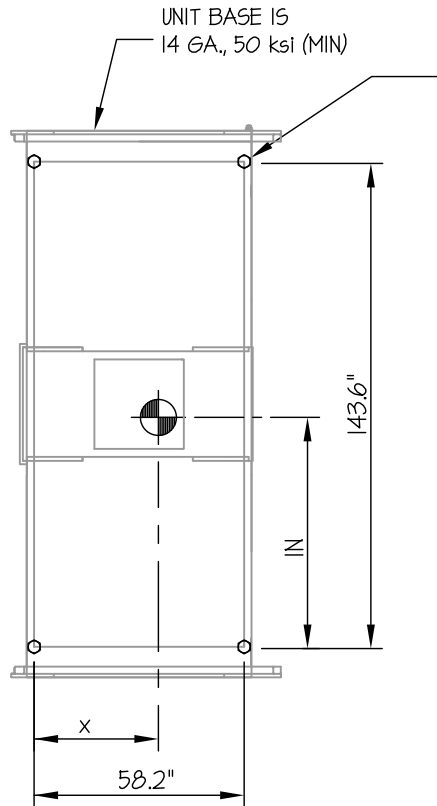

USE 4- 3/4" ϕ (A449)
THREADED RODS THRU
FLOOR (HOLES TO BE 1/16"
LARGER THAN BOLT
DIAMETER MAX.) TO STRUT
OR OTHER SUPPORTING
MEMBER (BY ENGINEER OF
RECORD)

MODEL	WEIGHT (LBS)	DIMENSIONS (INCHES)			BOLT FORCES	
		\bar{X}	\bar{Z}	\bar{Y}	T_u (LBS)	V_u (LBS)
*CS750-L	6700	33.4	63.1	53.5	4857	2884
CS750-M	5800	34.2	69.6	55.3	4058	2556
CS750-S	5400	32.6	72.4	59.8	4042	2269

* THIS MODEL REPRESENTED IN CALCULATION BELOW

PLAN AT BASE

$T_u = 4857$ LB/BOLT (MAX)
 $V_u = 2884$ LB/BOLT (MAX)



A P P R O V E D
Fixed Equipment Anchorage
Office of Statewide Health Planning and Development



OPA-2750-10

Pre-approval Program Manager:
Anthony R. Pike
(916) 440-8470

P E N D I N G

Reviewed By: Anthony R. Pike

9/22/11

BELIMED, INC.

CS750 CART WASHER TANKS

DES. **J. ROBERSON**

JOB NO. **11-1170**

DATE **9/22/11**



SHEET

8

OF **8** SHEETS

SEISMIC ANCHORAGE

CONCRETE SLAB ON METAL DECK

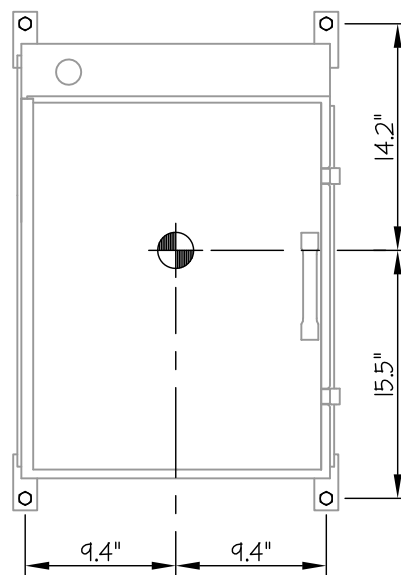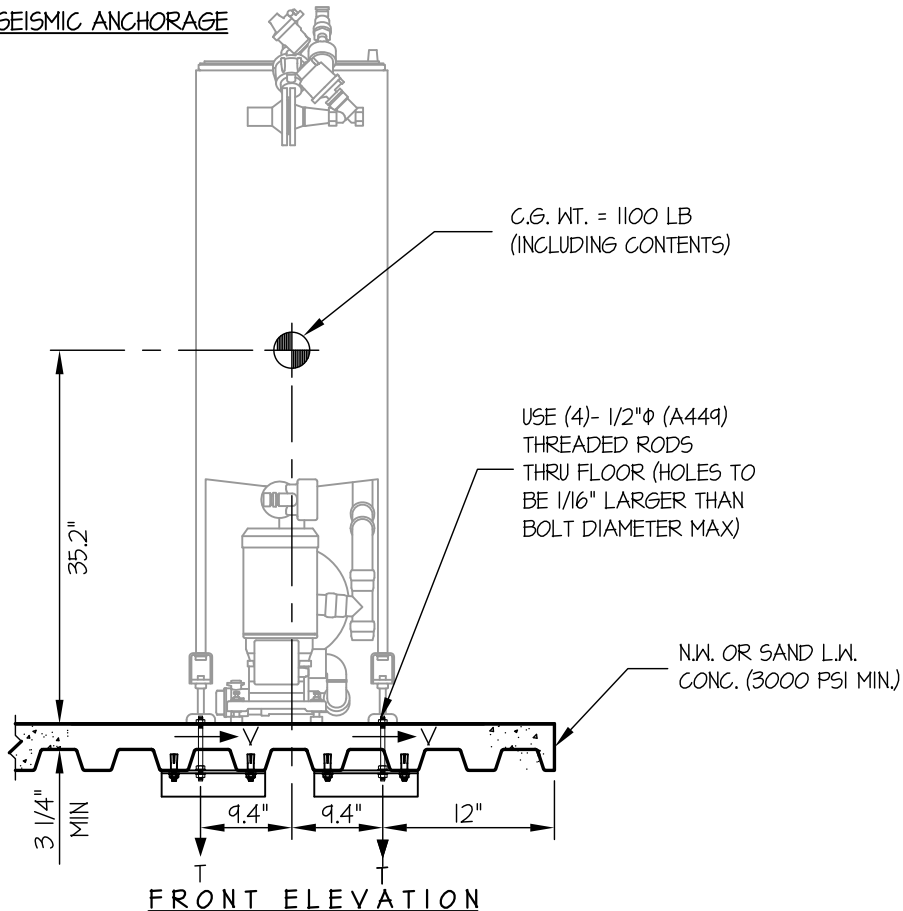

$T_u = 1762 \text{ LB/BOLT (MAX)}$
 $V_u = 431 \text{ LB/BOLT (MAX)}$



A P P R O V E D
Fixed Equipment Anchorage
Office of Statewide Health Planning and Development



OPA-2750-10

Pre-approval Program Manager:
Anthony R. Pike
(916) 440-8470

P E N D I N G

Reviewed By: Anthony R. Pike

9/22/11